



# vista®

for...

COCONUT WATER



## HunterLab

As a natural product grown in different regions of the world, the color of coconut water is impacted by the amount of pulp in the sample, and the color of that pulp. It is important to not only measure the color, but also the haze of the product. Haze% provides the best correlation to human perception of transmitted color. Vista® simultaneously measures visible-range transmittance color and haze% of coconut water in a single measurement. Vista® includes easy to navigate embedded EasyMatch® Essentials software and modern data communications. With its capacitive high-resolution color touchscreen display, Vista® sets a new standard for transmission color and appearance analysis of coconut water.



**Simultaneously measure the transmission color and haze% of coconut water**



**Haze% conforms to ASTM D1003 Procedure B, providing the best correlation to visual assessment**



**Operator friendly one-touch standardization – no accessories required**



**Embedded EasyMatch® Essentials software – no PC required**




**Smart user interface and data management**



measure  
**color**  
the way your  
**eye** sees it

## The Best Solution for Color Measurement of Coconut Water

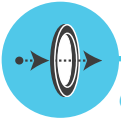

### Get instant color and haze results in a single measurement

Vista® offers this unique ability without the need to change software settings or change the sample container. Both color and haze% can be measured simultaneously in square cuvettes, rectangular cuvettes, and even round vials. This is as simple as it gets.



### Smart user interface and data management

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, and built-in email.



### Embedded EasyMatch® Essentials software – no PC required



The onboard EasyMatch® Essentials software is intuitive and easy to navigate. It offers applications to streamline production, allowing users to create their own workspaces to meet their specific application workflow needs. EasyMatch® Essentials comes preloaded with most standard transmittance scales and indices widely used in industry.

EasyMatch® Essentials




### Operator-friendly one-touch standardization

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.

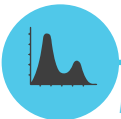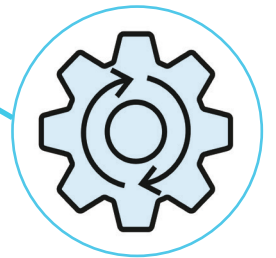

### Full spectrum, balanced LED illumination

Illuminating from 400-700nm, the balanced, full spectrum LED illumination is guaranteed to perform to original factory specifications for a minimum of 5 years of continuous measurements for years of worry-free operation

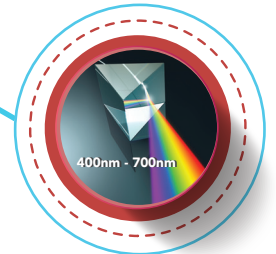

## Keys for the successful measurement of Coconut Water

Vista® offers a wide range of sample handling options accommodating round vials, ultra-micro and large format cuvettes, and even flow through cells. Contact our experienced technical staff today and we will help you to determine the right accessories to be successful.

