



Aeros®

Pat. #US 10,352,767 B2

for...

PHARMACEUTICAL PILLS



With new innovative SMART features like Auto Height Positioning, a rotating sample platform providing the largest sample area measurement in the world, large touch-screen display, and smart communications, the innovative Aeros® is the best way for you to measure the color of pharmaceutical pills of all kinds.

Aeros® combines patented non-contact measurement methodology with Auto Height Positioning, a rotating sample platform, and embedded EasyMatch® Essentials -our exclusive on-board color quality control software that offers all the color indices and metrics you need including CIE L*a*b*, CIE LCh, Hunter L,a,b and more.



Non-contact measurement eliminates and reduces sample cleanup time



Smart Sensor Technology with Automatic Height Positioning




Rotating sample platform provides the world's largest measurement area



Smart user interface and data management



Embedded EasyMatch® Essentials software – no PC required

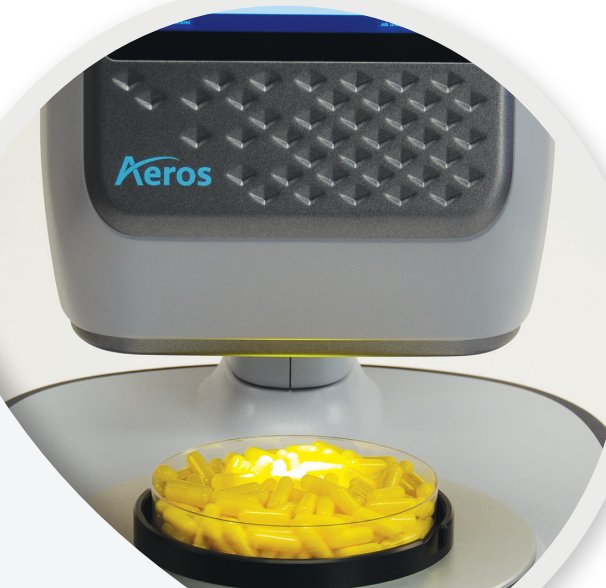


measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Pharmaceutical Pills



Non-contact color measurement

Aeros® measures pharmaceutical pills without touching them, simplifying sample presentation and eliminating possible contamination of either sensor or sample through contact. An added benefit is that sample clean-up time is significantly reduced.

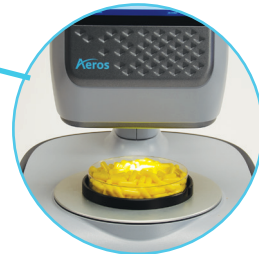


The world's only Smart Sensor Technology

The Smart Sensor Technology with Automatic Height Positioning simplifies operator involvement by determining the distance between the sensor and sample, and automatically adjusting for optimal measurement.



Dynamic rotating sample platform provides the world's largest measurement area

Simply place the sample and hit measure - samples automatically rotate beneath the smart sensor - 5 seconds, 35 measurements, 27.5 square inches of sample.

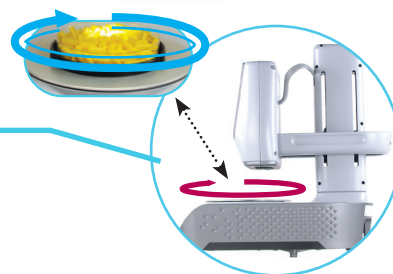

Smart user interface and data management

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, and built-in email.



Embedded EasyMatch® Essentials software – no PC required

The onboard EasyMatch® Essentials software is intuitive and easy to navigate. It offers applications to streamline production, allowing users to create their own workspaces to meet their specific application workflow needs. EasyMatch® Essentials comes preloaded with most standard reflectance scales and indices widely used in industry.

Keys for the successful measurement of Pharmaceutical Pills



Out of the box, Aeros® includes everything needed to be successful measuring pharmaceutical pills including two Sample Trays measuring 30.5 cm (12 in) Dia x 2.54 cm (1 in) Depth and 15.2 cm (6 in) Dia x 2.54 cm (1 in) Depth, along with Aeros® Quick Start Guide and USB flash drive.

